INSTRUCTION MANUAL

INDIVIDUEL PRO X

Aesthetic Professional
Home Device

Individuel
Genève

Dear Valuable Customer

Thank you for joining Individuel Genève

This is your smart move to your new beauty routine!

For consultation and assistance – contact us via:

- WhatsApp / phone: +41-76-673-6067
- Email: info@individuelcosmetics.com

Remember – REAL BEAUTY RUNS DEEP!

Individuel
Genève

Index			
Language	Page	Language	Page
English	3	Spanish	68
French	18	Korean	85
German	35	Hebrew	102
Arabic	52	Norwegian	119

Individuel Pro X

Table of Contents

PRECAUTIONS	3
WHO SHOULD NOT USE THE DEVICE?.....	3
CONTRAINDICATIONS.....	4
WARNINGS.....	4
WHAT IS INDIVIDUEL PRO X.....	5
INTENDED USE	7
POSSIBLE ADVERSE EFFECTS	7
DEVICE DESCRIPTION	8
User interface description	9
USER INSTRUCTIONS	10
What to expect when using Individuel Pro X	10
Preliminary Test	10
Prepare your skin	10
Device Set-Up	11
Device operation	12
Recommended treatment regimen	13
Post treatment facial care	14
CLEANING & STORAGE	15
TROUBLESHOOTING	15
CUSTOMER SERVICE	16
DEVICE SPECIFICATIONS.....	16
LABELING	17

PRECAUTIONS

- Do not use this appliance near bathtubs, showers or other vessels containing water as it may lead to electrocution risk.
- Do not use the Individuel Pro X in the presence of electronic monitoring equipment (e.g. cardiac monitors), as it may lead to their improper operation.
- Keep out of reach of children and pets.
- Children must not play with the appliance. Cleaning and maintenance by the user must not be carried out by children without supervision.
- Do not use the device if any part of the device is missing, damaged or loose, if the device is not working properly if you see or smell smoke.
- Before using the device, make sure the device is clean, dry and intact. Cleaning instructions are detailed in the "Cleaning & Storage" chapter.
- Do not use excessive pressure on the skin during treatment. Apply enough pressure for sensing the device's vibration according to your own comfort.
- Do not use Conductive Gel/cream/serum after its expiration date
- Use the device only on the allowed treatment areas and for the intended use.
- Use the device only according to the instruction for use and for its intended use.
- This appliance contains batteries that are non-replaceable.
- This appliance contains batteries that are only replaceable by skilled persons

WHO SHOULD NOT USE THE DEVICE?

Carefully read the information below before you begin a treatment with the device. If you still have any doubts or are unsure that the device is safe for your use, please consult with your physician.

CONTRAINDICATIONS

Do NOT use Individuel Pro X device if you:

- Are pregnant or attempting to become pregnant.
- Have any metal implant or an active implant such as a pacemaker, insulin pump or similar devices since it may cause interference.
- Keep out of reach of children and pets.
- Suffer from heart disease, such as arrhythmias.
- Have metal braces on your teeth
- Suffer from epilepsy or seizures.
- Suffer from blood coagulopathy or excessive bleeding or bruising.
- Have skin cancer, on areas at risk of cancer or in proximity to such areas.
- Have open sores or wounds in the treatment area.
- Tend to have rosacea or if you have problematic or sensitive skin.
- Are allergic to the Conductive Gel/cream/serum ingredients or sensitive to contact with chrome.

WARNINGS

Do NOT use the device:

- Over skin outbreaks, swollen or inflamed skin or lesions or if you have active skin disease (such as eczema, psoriasis and rash) in the treatment area.
- Over the neck and chest, since it may cause muscle spasms resulting in closure of your airway and heart rhythm disturbance.
- Over painful areas or areas that cannot sense touch. If you have painful or touch insensitive areas, consult with your physician before using the device.
- Without applying conductive gel/cream/serum/cream/serum before use.
- Individuel Pro X and the Conductive Gel/cream/serum are for external use only. Avoid direct contact with eyes, lips, mouth and nostrils. In case of contact, wash immediately with water.

Do not use Individuel Pro X while driving, sleeping, operating machinery or bathing/ showering.

- The appliance may be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge provided that they are under supervision or have been instructed in the safe use of the appliance and that they understand the hazards involved.
- If you have any uncertainty regarding any other medical condition not listed here, or if you have treated the skin with a prescribed drug or underwent any invasive procedure (such as injections or surgical procedures) within the past 3 months, consult your physician before using Individuel Pro X.
- Do not use the device while charging.

The Individuel Pro X device is designed for use on intact skin of all types unless prohibited in the above warnings and contraindications section.

WHAT IS INDIVIDUEL PRO X?

The Individuel Pro X is a facial toning device that is intended for facial stimulation and is indicated for over-the-counter cosmetic use.

The Individuel Pro X produces microcurrent. The microcurrent activity is accompanied by a mechanical vibration. This Individuel Pro X unique combination is especially designed to improve facial skin tone, stimulate the skin for improved appearance and relax facial muscles.

When using the Individuel Pro X you will need to use the Individuel Pro X Conductive Gel/cream/serum. The Individuel Pro X also includes a unique replaceable cartridge. This cartridge will allow you to release the desired amount of Conductive Gel/cream/serum effortlessly. The Conductive Gel/cream/serum is especially designed to help in proper conductivity and pleasant gliding of the Individuel Pro X across your skin.

The Individuel Pro X is designed for optimal facial contact for all people. The device should be used on intact skin of all types unless prohibited in the above warnings and contraindications section

How does it work?

The Individuel Pro X is an electronic non-invasive handheld device with a special and unique, easy-to-interact design. It is comfortable and simple to use, to facilitate several minutes of self-treatment on different facial areas.

The Individuel Pro X Conductive Gel/cream/serum is delivered to the selected treatment region with a simple click, the amount adjusted according to your comfort! This saves you the need of prior preparations. The Conductive Gel/cream/serum will enable proper skin conductivity and help the Individuel Pro X glide smoothly across the skin.

The microcurrent is delivered from the Treatment Surface when the device is in contact with the skin. This will automatically turn on the mechanical vibration.

The electrical microcurrent, accompanied with vibration, are used for aesthetic purposes, offering several benefits, including:

The electrical microcurrent, accompanied with vibration, are used for aesthetic purposes, offering several benefits, including:

- Improvement of facial contour and increased skin and muscle firmness
- Improvement of circulatory blood flow to the skin, thus facilitating waste removal and delivering fresh nutrients into the skin.
- Improvement in connective tissue collagen, thus supporting and structuring the skin.
- Increase in ATP (cellular energy molecule) synthesis, thus contributing to facial stimulation, a process that is slowed down by aging.

The microcurrent, together with the mechanical vibration, facilitates a healthier, more vibrant appearance!

 Please make sure to carefully read these Instructions for Use before you start your Individuel Pro X treatment. This guide, as well as an instructional video, can also be found at www.individuelprox.com

INTENDED USE

Individuel Pro X is intended for men and women who desire to improve their skin appearance, to tone their face and relax the muscles in the desired facial treatment areas.

Individuel Pro X is intended for use on the following facial areas: cheeks, around the eyes (periorbital region), forehead, around the mouth (perioral region), chin, neck, décolté and hands.

The Individuel Pro X device is designed for use on intact skin of all types unless prohibited in the above warnings and contraindications section.

Note: Individuel Pro X is especially designed for your own personal treatment! It is intended for a single user and home use only.

POSSIBLE ADVERSE EFFECTS

When using Individuel Pro X device according to the instructions for use, adverse effects are rare. However, every aesthetic procedure, including the usage of home-use devices, may cause some degree of adverse reactions (as listed in the table below). The possible adverse reactions are temporary, rare, and resolve spontaneously within 24 hours.

Severity of Adverse Effect	Frequency of Adverse Effect	Adverse Effects
Mild	Common	Local pain Skin irritation Skin redness Skin sensitivity Stinging Allergy or sensitivity to the gel/cream/serum
Moderate	Rare	Mild skin burns Excessive redness, rash and swelling Headache
Severe	Extremely rare	Severe burns Subcutaneous bleeding or bruising Infection Blistering Severe pain

In case adverse effects appear, stop using the Individuel Pro X device. Please consult your physician if necessary.

DEVICE DESCRIPTION

The Individual Pro X device set consists of the following components:

Hand-held device

Drop button: Conductive
Gel/cream/serum release
button

Cartridge level indicator

Remove Cartridge button

Conductive
Gel/cream/serum

Treatment Surface
(Large)

Battery
indicator

Electrical
outlet

On/Off
button


Hand Conductive
Surface

Additional device accessories

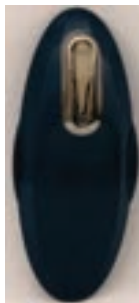1. Large Treatment
Surface

2. Small Treatment
Surface

User interface description

Button	Description
 <p>On/Off button</p>	<p>Press the button to turn On/Off the device. The device will turn off automatically when left unattended for 3 minutes.</p> <p>Green light – Device is on (microcurrent & vibration works when the Treatment Surface & the Hand Conductive Surface are in proper contact with the skin).</p> <p>Green light will appear:</p> <ol style="list-style-type: none"> 1. After turning on the device. 2. After pressing the Cartridge Removal button.  <p>When the Red light appears, do not remove cartridge status.  The light will turn back green after pressing the Cartridge Removal button</p>
 <p>Drop button: Conductive Gel/cream/serum release</p>	<p>Press  to apply the Conductive Gel/cream/serum on the treatment area.</p>
 <p>Cartridge Removal button</p>	<p>Press the Cartridge Removal button in order to remove the cartridge.</p>
<p>Cartridge level indicator</p>	<p>Indicates the level of remaining Conductive Gel/cream/serum in the cartridge:</p> <ul style="list-style-type: none"> ● Full cartridge (4 lines) ● ¾ cartridge (3 lines) ● ½ cartridge (2 lines) ● ¼ cartridge (1 line) ● empty (no light appears)
<p>Hand Conductive Surface</p>	<p>A secondary electrode that closes the microcurrent circuit with your palm. Place the palm of your hand on the Hand Conductive Surface to start treatment.</p> <p>It is recommended to use a moist hand for improved conductivity.</p>
<p>USB cable</p>	<p>USB cable is used to recharge the battery of the Individuel Pro X.</p>
<p>Battery Indicator</p>	<ul style="list-style-type: none"> ● Green light indicates that the battery is charged. ● Red light indicates that the battery is low and needs to be recharged. ● When charging is complete the light will turn off. ● The battery lasts for approximately 2.5 hours of constant work.

USER INSTRUCTIONS

What to expect when using Individuel Pro X

The use of Individuel Pro X device for facial stimulation for aesthetic purposes is simple and quick.

During Individuel Pro X treatment, it is normal to experience:

- Sensations of warmth and vibration, created by the microcurrent energy and the vibration motor.
- Some mild red or pink color on the skin and irritation sensation, during and just after using Individuel Pro X
- Motor noise during movement of the cartridge into/out of position.
- Vibration motor noise during treatment.

Preliminary Test

Perform a preliminary test to make sure your skin's reaction to the Individuel Pro X device is normal. Perform a preliminary test on a small patch of skin in an allowed facial area, prior to the first treatment with the Individuel Pro X device.

1. Prepare your skin (as described in the "Step 1: Prepare your skin" section) and set up the Individuel Pro X as indicated in "Step 2: Device set up".
2. Operate the device for 1 minute as detailed below in "Step 3: Device Operation".
3. If you feel uncomfortable or see adverse reactions (described in the "Adverse Reactions" section), stop the test immediately.

Wait 24 hours and examine the tested area. If the skin in the tested area shows no response or only mild redness, you can proceed with a full treatment. If you see extreme redness, swelling or blistering, do not proceed with the treatment.

Step 1: Prepare your skin

The facial skin is a sensitive area that requires special care and caution. Make sure you follow all instructions and safety precautions carefully, in order to avoid adverse reactions.

1. Make sure your skin is clean and free of makeup and jewelry. Oil residues or dirt might reduce the treatment efficacy.
2. Perform a preliminary test on a small patch of skin in an allowed treatment area (as described in the "Preliminary Test" section).

Step 2: Device Set-Up

What to expect when using Individuel Pro X

1. Remove all components from the box (hand-held device, treatment surfaces, conductive gel/cream/serum cartridge and USB cable).
2. Charge the Individuel Pro X for 6 hours prior to using it for the first time.
3. Remove the Conductive Gel/cream/serum Cartridge from its packing (see Figure 1)
4. Choose the required Treatment Surface (large or small):

Small Treatment Surface – around the eyes, around the mouth and chin.

Large Treatment Surface – full face treatment (including cheeks) neck, décolte and hands.

5. Place the cartridge inside the Treatment Surface with the cartridge tip facing the Treatment Surface pore (see Figure 2 & 3).
6. Insert the Treatment Surface loaded with the Conductive Gel/cream/serum Cartridge into the device and press firmly to snap it into place (Figure 4).

Figure 1 & 2 – Insert the Conductive Gel/cream/serum Cartridge tip inside the selected Treatment Surface pore.

Figure 3-Insert the Treatment Surface with cartridge inside the Individuel Pro X

Figure 4-Your Individuel Pro X is ready for use!

7. Before beginning the treatment, check your device: make sure the device and the Treatment Surface are clean, dry and intact
8. **Note:** If you see sharp edges or defects on the Treatment Surface, do not use the Treatment Surface and contact the Individuel Pro X Support Service.
9. The Individuel Pro X is now ready for use.

Step 3: Device operation

1. Turn on the device by pressing the On/Off button.
 - i. The On/Off button will turn light green.
 - ii. The cartridge level indicator will light-up
 - iii. The battery indicator will light green, meaning that the battery is charged.
2. Press the Drop button  to set the cartridge position:
 - i. The device motor will move until reaching the cartridge. The device motor will set itself automatically according to the amount of Conductive Gel/cream/serum left in the cartridge.
 - ii. The On/Off button  color will change from green to red.
 - iii. The cartridge level indicator will show 1-4 green lights, depending on the amount of Conductive Gel/cream/serum left inside the cartridge (See the “User Interface Description” section). When the cartridge is empty the indicator will not light up.
3. Place the device on the desired treatment area and use the Drop button once or several times, to apply the required amount of Conductive Gel/cream/serum until the desired treatment area is covered. Each press releases one drop.
4. Make sure you have applied the required amount of Conductive Gel/cream/serum. Keep your treated area moist with Conductive Gel/cream/serum throughout the entire treatment session
5. The palm of your hand must always touch the Hand Conductive Surface. Moisten your hand in order to enhance conductivity.

Note: If you can feel a gentle vibration, the microcurrent is being delivered.
6. Repeat points 3-5 until you have covered the entire treatment area of your choice.

Note: In case the device identifies insufficient contact with the palm of your hand, the current and vibration will be disabled.
7. Turn off the device by pressing the On/Off button or the device will shut down automatically after 3 minutes of non-use.

8. For Individual Pro X device maintenance and prolonged use, cleaning is required (As described in the "Cleaning & Storage" section).

9. To replace between Treatment Surfaces or to remove the Conductive Gel/cream/serum Cartridge, make sure that the On/Off button light is green.

i. If the On/Off button light is red, press the Remove Cartridge button  and wait until it becomes green

ii. Remove the cartridge from the Treatment Surface.

iii. Insert the Conductive Gel/cream/serum Cartridge to another Treatment Surface or replace the Conductive Gel/cream/serum Cartridge with a new one (according to "Step 2: Device Set-Up").

iv. Insert the Treatment Surface with the Conductive Gel/cream/serum Cartridge into the device (according to "Step 2: Device Set-Up").

10. If the Battery Indicator is red, the battery needs to be recharged. Plug the USB cable into the device and plug the other end of the cable into the USB source. While the device is charging, the indicator light will be green. When charging is complete the light will turn off.

Note: A full cartridge contains 15ml of Conductive Gel/cream/serum, which is enough for about 5-6 full facial treatments.

Recommended treatment regimen

Always perform a preliminary test prior to the first treatment in each facial area according to the "Preliminary Test" section.

Treatment Regimen	Treatment Protocol	Treatment Area and position	Treatment Surface
Full Face including Cheeks, neck, decolte and hands	<ul style="list-style-type: none">• Treat entire face for up to 10 min as shown in the picture.• Repeat daily according to treatment protocol		

Forehead	<ul style="list-style-type: none"> • Treat for up to 5 min as shown in the picture. • Repeat daily according to treatment protocol 		
Around the Eyes	<ul style="list-style-type: none"> • Treat for up to 5 min as shown in the picture. • Repeat daily according to treatment protocol 		
Around the Mouth and Chin	<ul style="list-style-type: none"> • Treat for up to 5 min as shown in the picture. • Repeat daily according to treatment protocol 		

Note: A few (2-10) drops may be required for each treatment area.

Note: Make sure the device is always in contact with the skin. When there is no skin contact, you will not feel the mechanical vibration and microcurrent will not be delivered.

Post treatment facial care

Following treatment, avoid the use of abrasive or harsh products for two days.

- In cases where superficial crusting or scabbing occurs, do not rub, scratch or remove the scab. Cleanse the area gently and pat dry.
- Make-up may be applied immediately after the treatment.

CLEANING & STORAGE

- Treatment Surfaces should be cleaned after each treatment:
 - If the On/Off button  is red, press the  Remove Cartridge button. When the On/Off button  turns green, the Treatment Surface can be removed.
 - Remove the Conductive Gel/cream/serum Cartridge from the Treatment Surface.
 - Store the Conductive Gel/cream/serum Cartridge until the next treatment in its package in a dry place.
 - Rinse the Treatment Surface thoroughly.
 - Make sure the Treatment Surface is dry (you may use cloth or paper towel to wipe the water) and store it.
- Make sure the device is turned off and unplugged, before cleaning it of residual Conductive Gel/cream/serum. Use a soft, damp cloth or paper towel. Make sure that all areas including the Treatment Surface and buttons are thoroughly cleaned.
- Make sure the Individuel Pro X device is completely dry before storing it; never rinse the device with water as it may lead to electrocution risk.
- Store the Individuel Pro X in a closed and dry place. It is recommended to keep it in its original box.

TROUBLESHOOTING

If you experience any problems with your Individuel Pro X device please refer to the information below:

If the Individuel Pro X device does not turn on	<ul style="list-style-type: none">• The battery might be discharged. Try to charge the device with the USB cable.
If the Individuel Pro X does not deliver current (no mechanical vibration):	<ul style="list-style-type: none">• Make sure that the palm of your hand is in full contact with the Hand Conductive Surface (a slightly moisten palm is recommended).• Make sure that the treatment surface is in contact with the treated skin that has conductive gel/cream/serum on it.• If your palm is not properly applied and/or the treatment surface is not in contact with your skin, the current will not be delivered.
General error:	<ul style="list-style-type: none">• If there is dirt or debris stuck to the electrode, you should clean them in accordance with the “Cleaning & Storage” section above. If you cannot• Reset the device – turn it off or wait until it shuts down, wait 10 sec and then turn it on again.• If you still experience problems with your device, please contact our customer service.

CUSTOMER SERVICE

For customer service, please see contact details on the inner front cover and refer to our website: www.individuelprox.com.

DEVICE SPECIFICATIONS

The Individuel Pro X device is powered by (applied with continuous current) $5\pm 1[V]$.

Applied current: 5-100 [μA]

Mechanical vibration- Implemented by treatment surfaces 3.4 [G]

Chrome plated Treatment Surface:

- Small Treatment Surface: ~2 [cm²]
- Large Treatment Surface: 22 [cm²]

USB Cable

Device Size and Weight: 192.6X48.9X64 [mm] (LXWXH), 120 [gr]

Storage Instructions:

Storage and transport conditions	
Temperature*	-10°C to 70°C (14°F to 158°F)
Humidity	10% to 90% RH at 55°C/131°F
Pressure	Non-condensing
Operational conditions	
Temperature*	10°C to 40°C (50°F to 104°F)
Humidity	30% to 70% RH
Pressure	3,000m.

Use Individuel Pro X only with the supplied external power and accessories. For your safety, never use DC adapters which are not approved by the manufacturer.

* In case of storage or transportation at extreme conditions, allow the Individuel Pro X device 30 minutes for thermal stabilization

LABELING

Model: DEP100XXYYZZZ

Power: 5 [V], 0.25[A]

Chargeable Li-Ion Battery: 3.7V 190mAh

	Class III insulation
	Type BF equipment
	Consult Instructions for use
	Waste Electronic and Electrical Equipment
	Our products and packaging are recyclable, do not throw them away!
	Keep away from eyes
	Please read all warnings and safety precautions prior to use
	Do not use the Individuel Pro X device on sensitive tanned skin or after sun exposure
	Do not rinse Individuel Pro X under running water, as it may cause device malfunction. Keep dry.
IP22	IP2X: Object of 12.5 mm diameter, shall not penetrate. IPX2: Vertically falling drops shall have no harmful effects.
CE	
UK CA	

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